

HILL 90C

Manual and Air Drop Tables




Hill Laboratories
COMPANY

HILL 90C

Manual Drop Table

Designed with both you and your patient in mind, the standard 90C electrically adjusts in height from 22" to 35" providing the ideal working height for you and easy access for your patients. The 90C is also versatile—choose from the many available options to customize the 90C for your specific needs.



Pictured with several options

Shown in Suede Upholstery with Beige Base.

▼ STANDARD FEATURES

- Armrests, 11" face cut-out and paper roll
- Electrically-controlled height - 22" - 35"; Optional 20" - 30"
- Width - 24"; Length - 5'10" with drops; 6' without drops
- Foam Top - 2 1/2" Medium (standard); Soft or Firm
- Choice of 20 Colors for the top cushions
- Beige, Gray or Black Base, enclosed for safety
- One-year warranty on all parts and labor



The 90C has electric adjustable height from 22" to 35" and available in three base colors; Beige, Gray and Black.

TOP PICK

The model shown is one of our most popular configurations with optional tilt, raised and dual-drop headpiece, thoracic drop, lumbar drop, pelvic drop and slide-out ankle support. For even greater speed and convenience, consider the Air90C.

AIR90C

Air-Drop Table

Crisp, smooth and impressively fast.

Equipped with the Hill Air-Drop System*, the Air90C allows doctors to select one or up to four drops simultaneously from the convenient, side-mounted, control panel. Cushions are activated instantly by tapping one of the foot strips.

Shown in Cinder Upholstery with Black Base and added options.

** The Hill Air-Drop System, with air-powered features, is only available on the Air90C.*

*Dual Power Foot Strips
Reset the Drops with a tap!*

STANDARD FEATURES

- Height - 21 1/2" to 31" (higher starting heights available)
- Rocker Foot Pedal and 2 Power Foot Strips
- Width 24", Length 6'3"
- Armrests, 13" Face cut-out, Paper Roll, Cutter
- Tilting Headpiece, 30° positive and negative tilt
- Air Cervical Dual-Drop Headpiece
- Air-Powered Thoracic Spring Breakaway and Release
- Air Thoracic, Lumbar and Pelvic Drops
- Drop-Away Foot Section
- Choice of 20 Colors for the top cushions
- Beige, Gray or Black Base, enclosed for safety
- One-year warranty on all parts and labor

Select one or up to four drops simultaneously from the Control Panel. Optional dual-sided control box available



90C OPTIONS

Made-to-Order

▶ MANUAL OR AIR DROPS


The standard 90C can be equipped with up to 4 manual drops including cervical (straight-down and forward-motion), thoracic, lumbar and pelvic. The Air90C comes standard with all 4 air drops.

▶ THORACIC SPRING BREAKAWAY


The Thoracic Spring Breakaway (with release) allows for a thoracic recoil movement, features adjustable tension and locks at any point down to 5 1/2"; with drops - 4 1/4". This feature cannot be combined with a thoracic or lumbar drop. The Air90C comes standard with an air-controlled breakaway and pressure pedal.

▶ HEADPIECES

Add functionality to your headpiece by combining several headpiece options including Tilt, Dual-Drop, Raised, Flexion, Axial-Motion, Lateral and Translation (CBP). The Air90C comes standard with a Tilt and an Air Dual-Drop combination headpiece.

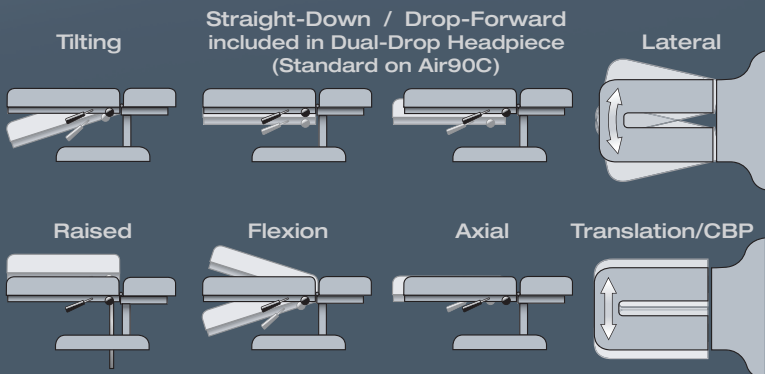

▶ MANUAL CERVICAL FLEXION

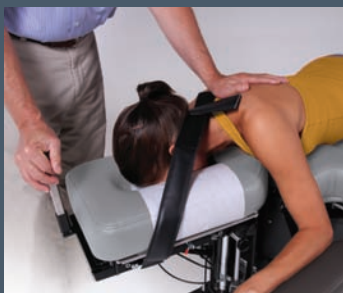


The Flex-Headpiece comes with the raised headpiece, occipital strap and removable tiller bar. The headpiece unlocks and a gas-spring gently compresses and recoils as you control the flexion, lateral and long-axis movement.

All of the options below can be added to the 90C. Most are already included on the Air90C as standard features. These are indicated by the blue "Air" symbol. **AIR**


See Air90C page for a complete list of features.

Thoracic/Lumbar Options

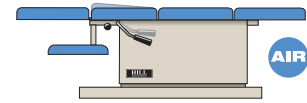


Thoracic Drop

The 90C manual drop is activated with an easy-lift handle. The Air90C is activated with Foot Strips.

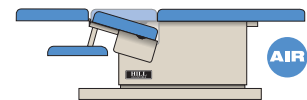


Thoracic Release

Hinges and locks at any angle down to 5 1/2 inches.

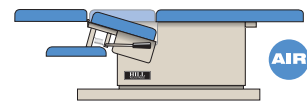

Thoracic Breakaway

Allows for a thoracic recoil movement and features adjustable tension. The Air90C comes standard with an air-controlled breakaway and pressure pedal.

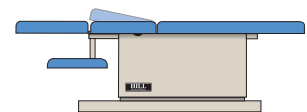

Raised Thoracic

Tilts and locks up to 6 inches. (not available with spring breakaway or Air90C)



Lumbar Drop

The 90C manual drop is activated with an easy-lift handle. The Air90C is activated with Foot Strips.

any are

atures.




Ask about our flexible lease-to-own options available for 1, 2 and 3 year financing. Your 90C may also be eligible for a tax credit totaling nearly 1/2 the purchase price.

Visit HillLabs.com for details.

Pelvic Options


Other Options

Special Configurations

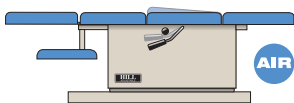


Pelvic Drop

The 90C manual drop is activated with an easy-lift handle. The Air90C is activated with Foot Strips.

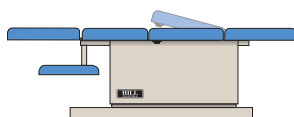

Raised Pelvic

Tilts and locks up to 5 inches.

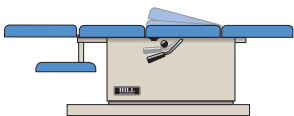

Raised Pelvic/Drop Combo

Combines the Raised Pelvic feature with a Pelvic Drop.

Standard pelvic cushion is notched; solid cushion upon request. All 90C drops are activated with an easy-lift handle on both sides of the table. All Air90C drops are selected using the control box and activated using the Foot Strips. Opposite-side control box is available.



Drop-Away Foot

9" foot section folds down for treating shorter patients. Standard on Air90C.




Slide-Out Ankle Support

Adjusts up to 8" in height and extends out 12" from the end of the table.




Adjustable Armrests

Elevating armrests for supine shoulder support – often used for massage.

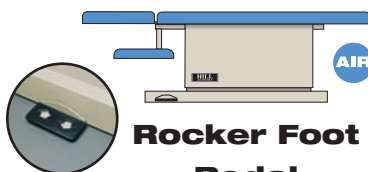

Rocker Foot Pedal

The Rocker Foot Pedal attaches to the base on both sides of the table to control the height. Standard on the Air90C.

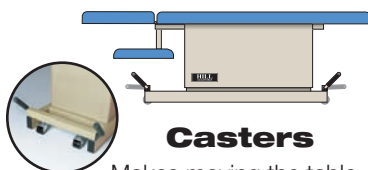


Castors

Makes moving the table easy. Engaged with a foot lever.



Extended Face Cut-Out

Extended face cut-out for treating children and shorter patients. Not available on Air90C. Other custom face cut-outs available.

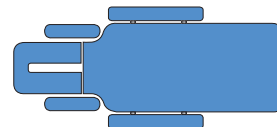

Removable

Width Extensions

3" x 21". Adds 6" to the width of the table. Used for massage therapy; often combined with adjustable armrests.



Activator*-Style Top

36" face cut-out and wedged foot section. 27" width, optional 30". Not available on Air90C.



Activator*-Style Top

Side view showing the wedged foot section - Model HA90AM.

The Activator-Style top is also available as a Hi-Lo model with soft start-up and landing plus automatic fold-down standing board.*

SPECIFICATIONS

- **Powered by an electric actuator.** 110 volts, 60 Hz, 5.0 amps.;
Air90C - 2.5 amps. 220 Hz - 50 cycle available at additional charge
- **Lifting capacity** - 450 lbs.
- **Travel time** from lowest to highest height - 15 seconds
- **Choice of colors** for upholstered top (see below)
- **Beige, gray or black telescoping base skirting** constructed of fire-rated acrylic thermo-plastic
- **Powder-coated** internal and external parts
- **Medium density foam top** - 2 1/2"; firm and soft density available
- **Standard width** - 24"; optional widths, 22", 27", or 30"
- **Range of height**
90C from 22" to 35"
- with thoracic release or breakaway option, 24" to 35"
- custom low height available
Air90C from 21 1/2" to 31"
- 23 1/2" to 33" available. Special heights upon request
- **Table Top length**
90C - 5'10"; with drop-away foot becomes 5'1"
Air90C - 6'1"; with drop-away foot becomes 5'4"
Slide-out foot section: removed; 5'1", inserted 5'4", extended 6'4"
- **Shipping weight** - approximately 250 - 300 lbs.
- **Warranty** - one year on all parts and labor

*The option shown for use with the ACTIVATOR technique is designed to serve the doctor who uses the ACTIVATOR® instrument or applies the "ACTIVATOR METHODS®" technique. ACTIVATOR and ACTIVATOR METHODS are trademarks of Activator Methods, Inc. and are used for identification purposes only. Equipment made or distributed by Hill Laboratories is not made or approved by Activator Methods, Inc.


VERY AFFORDABLE

Finance your Hill Laboratories table through our convenient and affordable lease/purchase plan. Visit HillLabs.com to calculate leasing options or call to speak personally with one of our representatives who will be happy to discuss payment options.

AVAILABLE COLORS



Below - Finesse® Colors. Upgraded material; foam-backed and puncture-resistant.

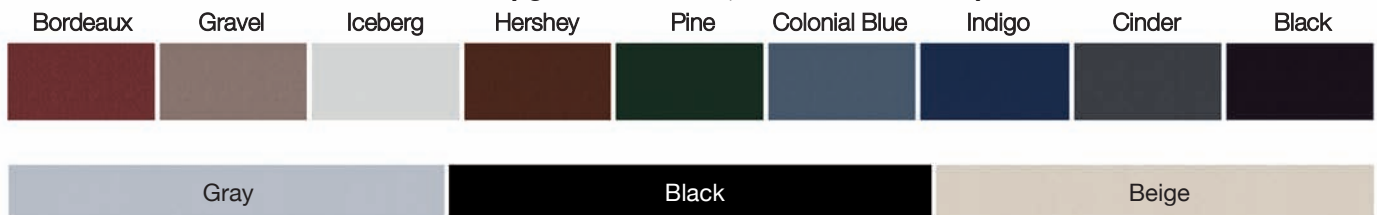

Upholster your table in durable Boltflex® with Prefix™, scientifically formulated to give superior protection from staining and abrasions. Finesse and custom colors are available at additional cost. Shades may vary from actual material. Color choices may change without notice.



P.O. Box 2028, 3 N. Bacton Hill Road, Frazer, PA 19355 • Toll-Free 1-877-445-5020
610-644-2867 • fax 610-647-6297 • www.HillLabs.com